



BICEPS Overview



Revision 1,
2018-10-02

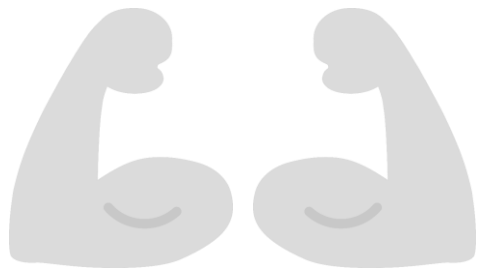

Orientation

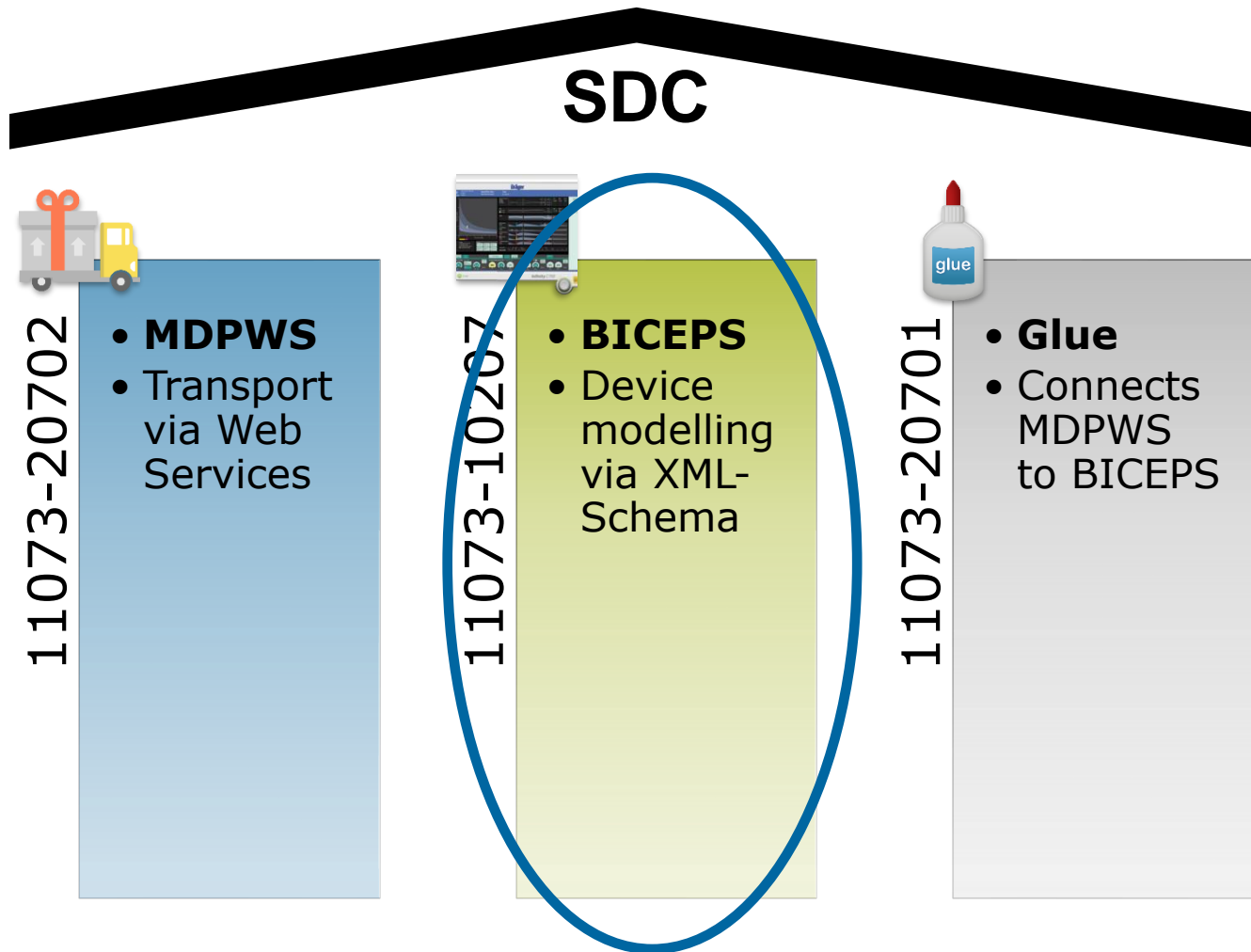


At a glance

- Official title: 11073-10207 – Domain Information & Service Model for Service-Oriented Point-of-Care Medical Device Communication
- Non-normative title: BICEPS = Basic Integrated Clinical Environment Protocol Specification
- Conceptual model based on ideas of 11073 classic DIM and SOMDA
- Semantic description of medical device capabilities and state information
- Compatible with ICE architecture led by MD PnP


BICEPS does not define means to convey data over a physical layer!

SOMDA


SOMDA = Service-Oriented Medical Device Architecture

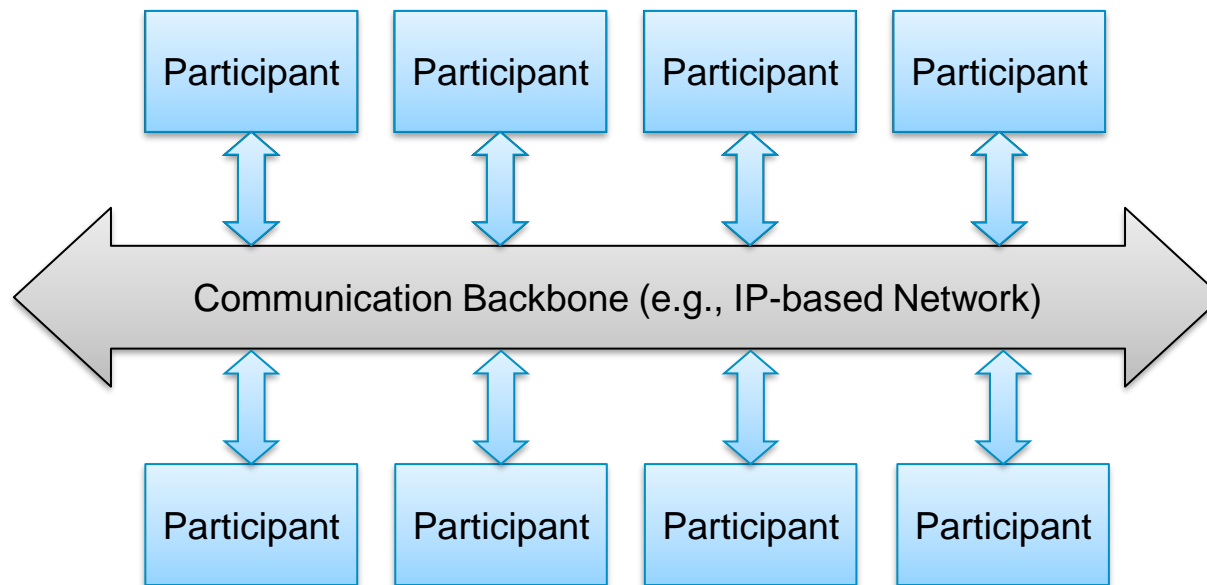

Wording

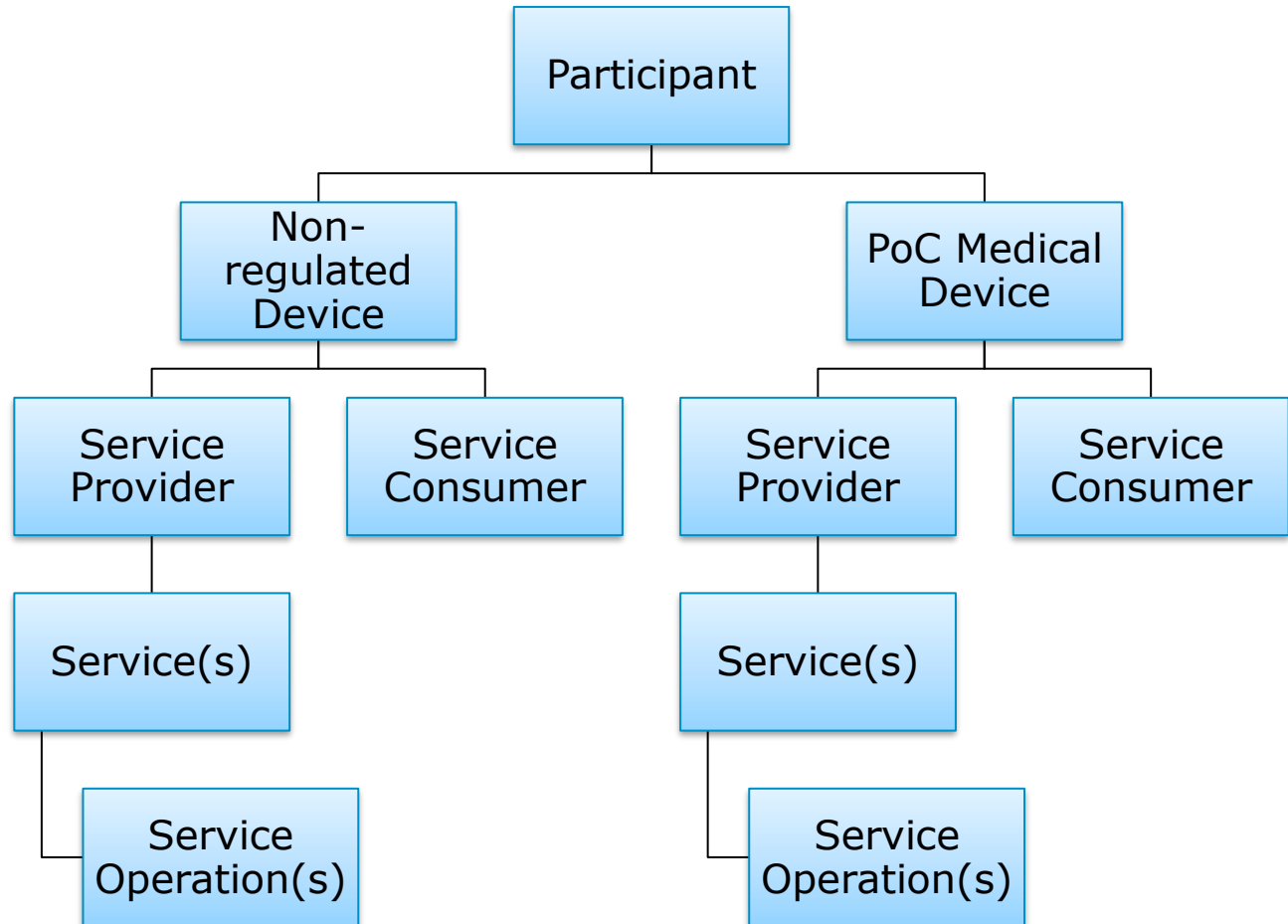


Component view

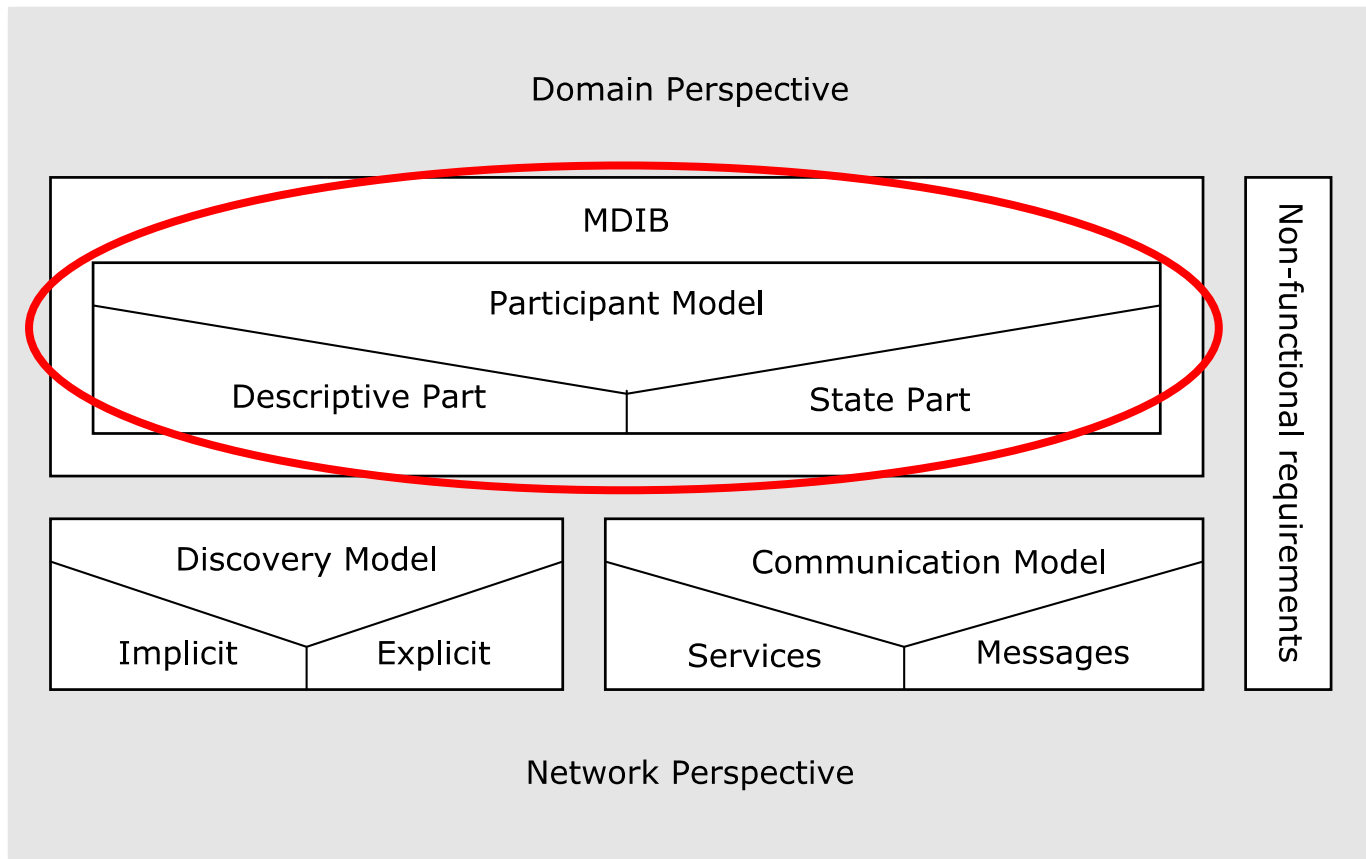


Participant model

MDIB

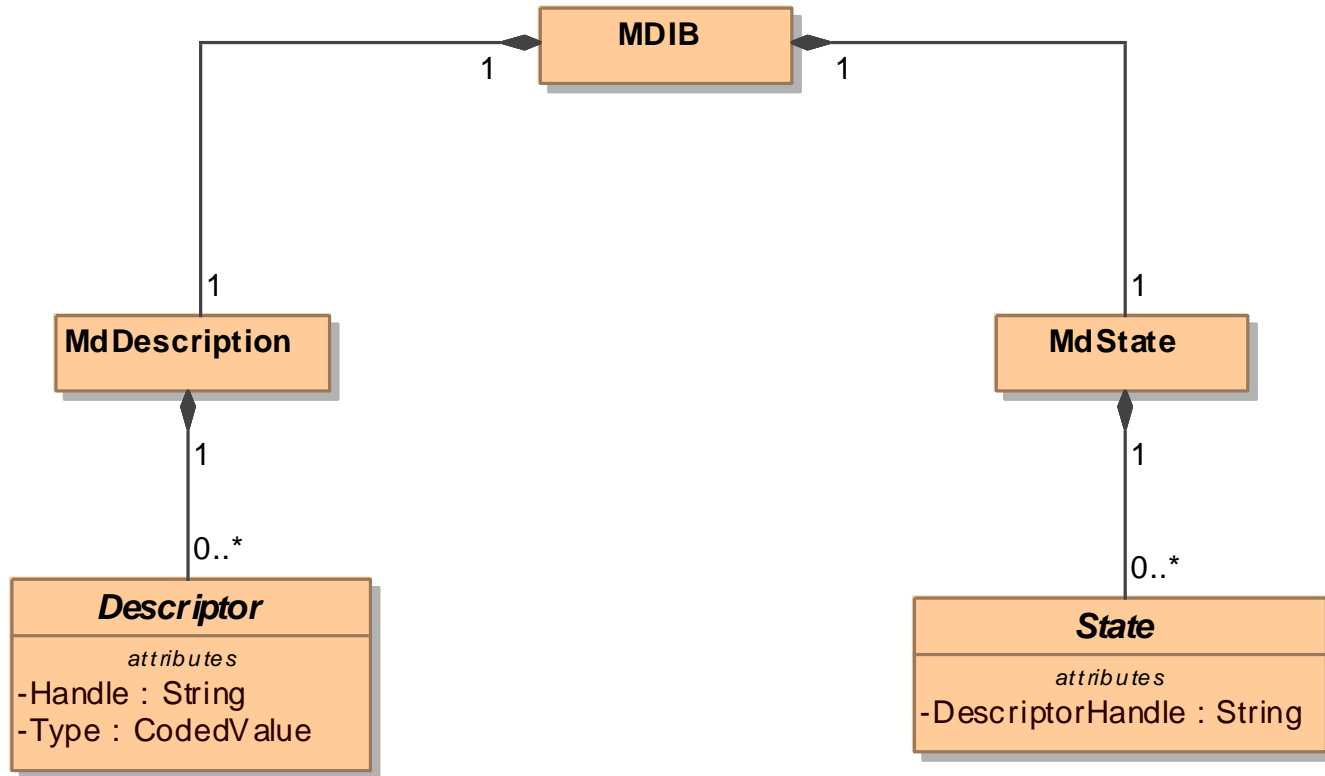


Participant model

Descriptive part

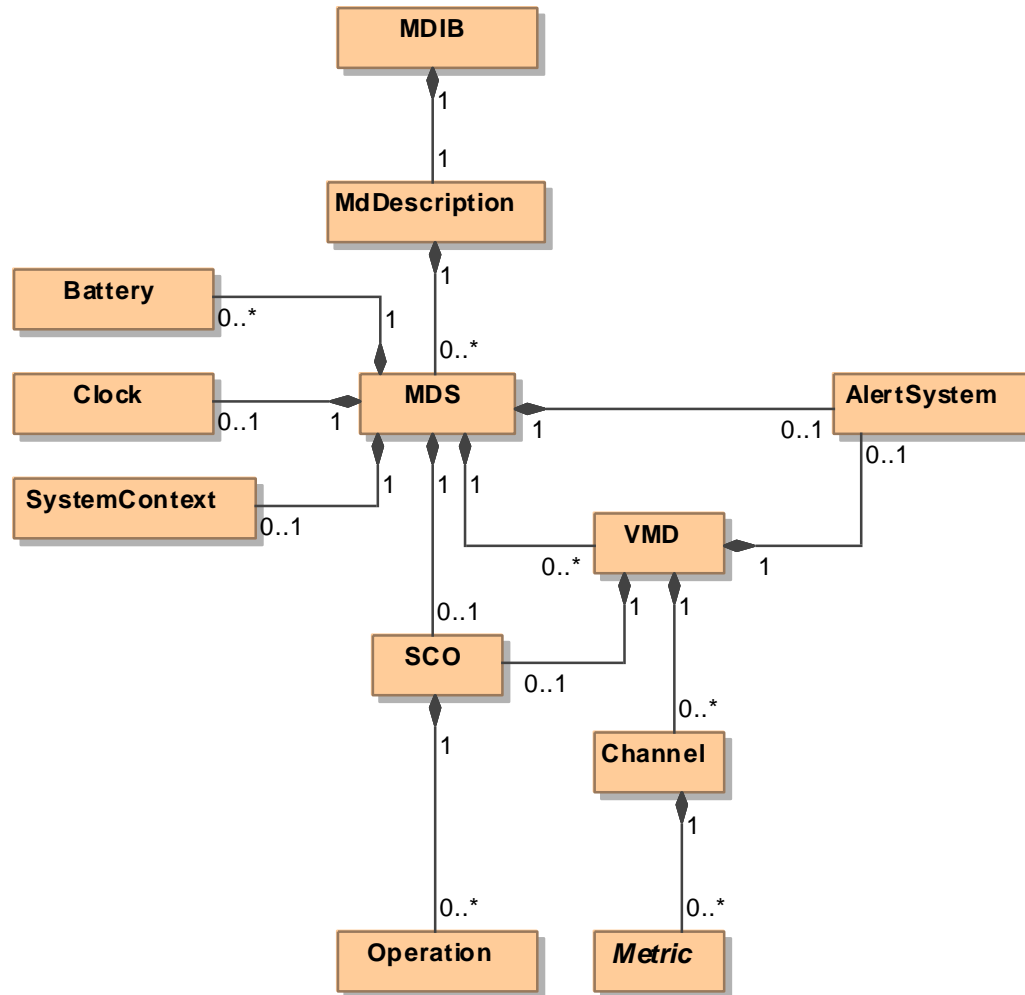

Participant model


Containment tree

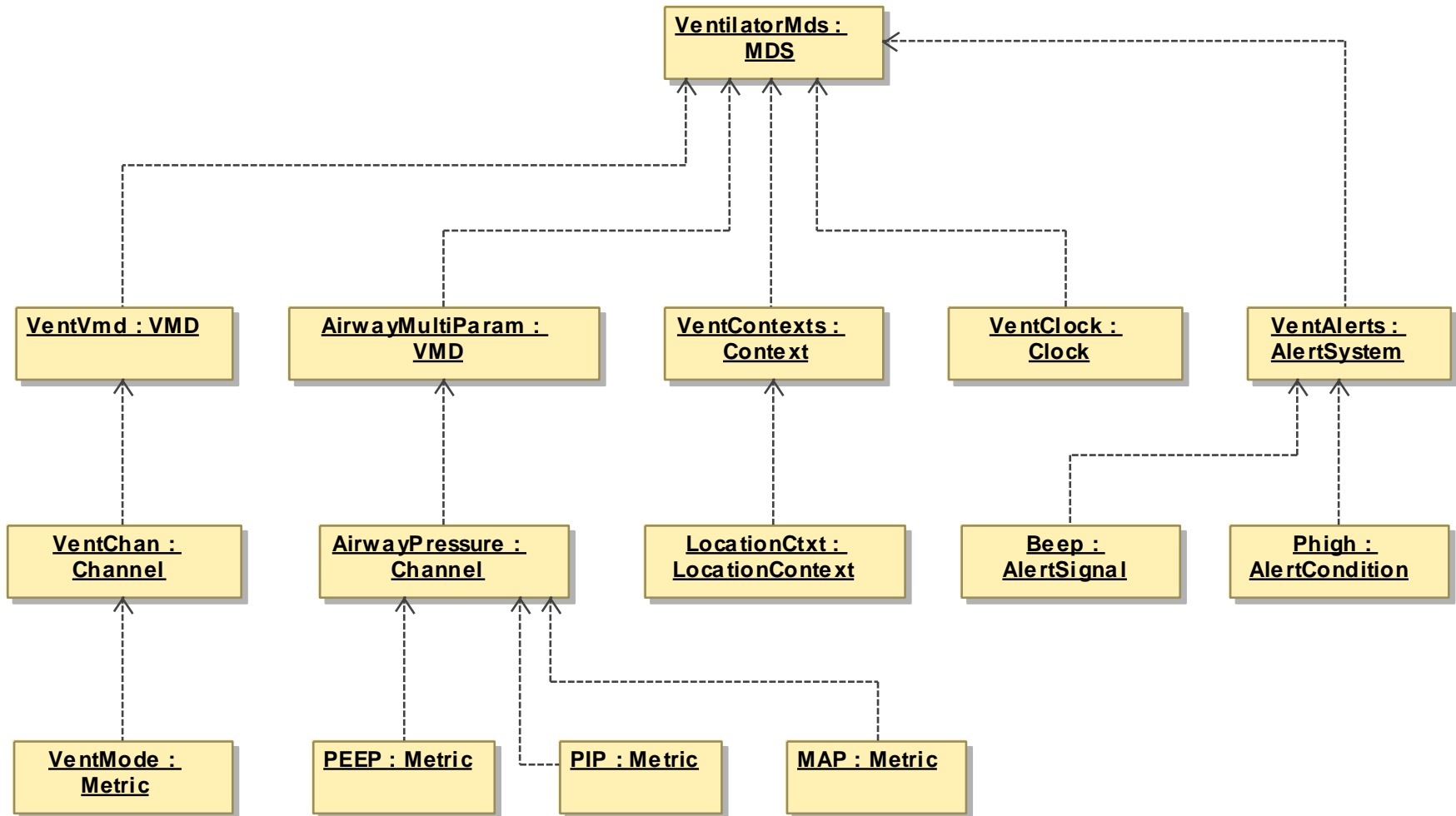


Participant model

State part

- Volatile state (of measurements, settings, contextual info)
- Might change frequently
- Single state: w/o handle; is identifiable by the descriptor handle
- Multi state: w/ handle; because not uniquely identifiable by the descriptor

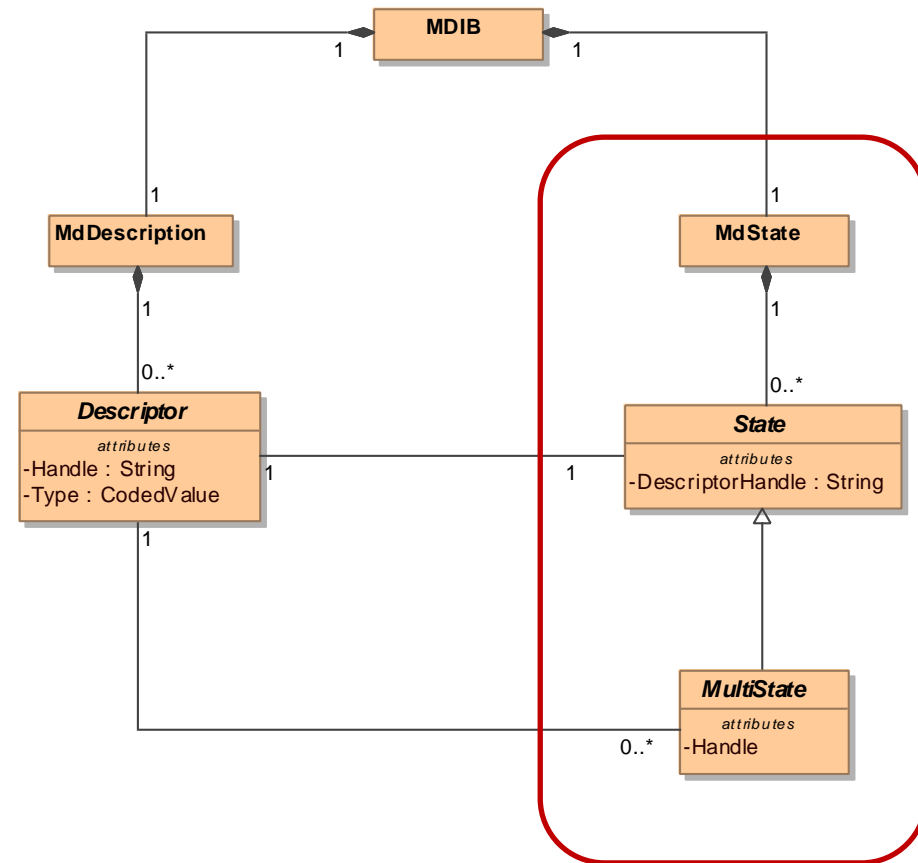

Component view

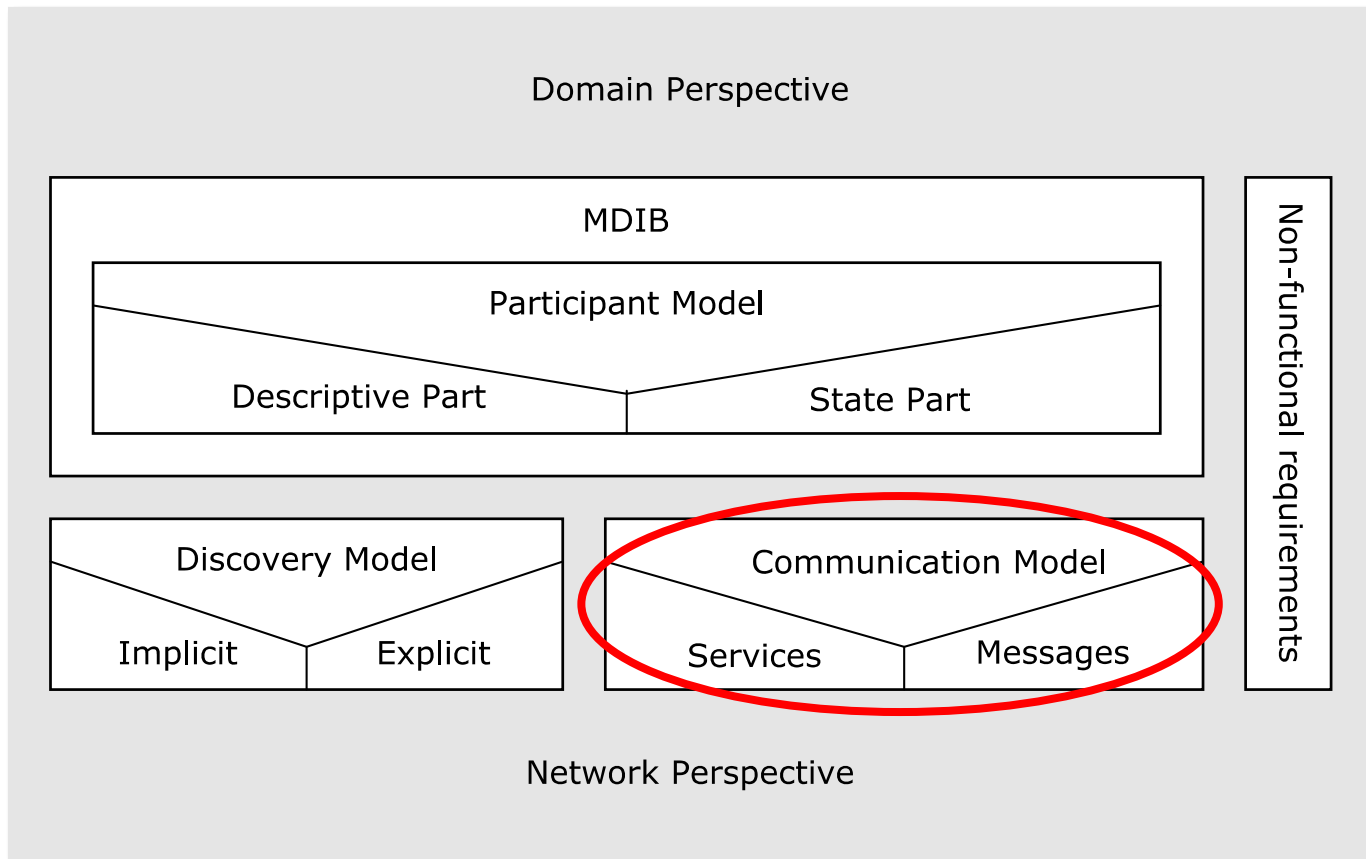

Communication Model

- BICEPS is based on Service-oriented Architecture (SOA)
 - In SOA, service operations are invoked to exchange messages
- BICEPS Communication Model defines
 - Set of services
 - Set of messages

Services


What is a Service?

- A service is an abstract construct in a Service-Oriented Architecture (SOA)
- A service is hosted at a service provider
- Service consumers invoke service operations of the services a service provider exposes to the SOA communication backbone
- In BICEPS, there are two kinds of service operations
 - Request-response
 - Driven by the service consumer
 - Notification
 - Event-driven by the service provider
- BICEPS is based on SOA principles and does not define an implementation
 - the Glue spec is a set of rules that ties together MDPWS and BICEPS

Services

Means to access the MDIB


- **GET**
mandatory, request-response of whole MDIB, or descriptive and state part separately
- **SET**
MDIB SCO for remote control
- **STATE EVENT**
State change events
- **DESCRIPTION EVENT**
Descriptor change events
- **CONTEXT**
Retrieve and set context states
- **WAVEFORM**
Retrieving streaming data
- **CONTAINMENT TREE**
Traverse containment tree and retrieve descriptors in finer granularity
- **ARCHIVE**
Access historical data
- **LOCALIZATION**
Retrieve human-readable descriptions if not stored in the MDIB directly

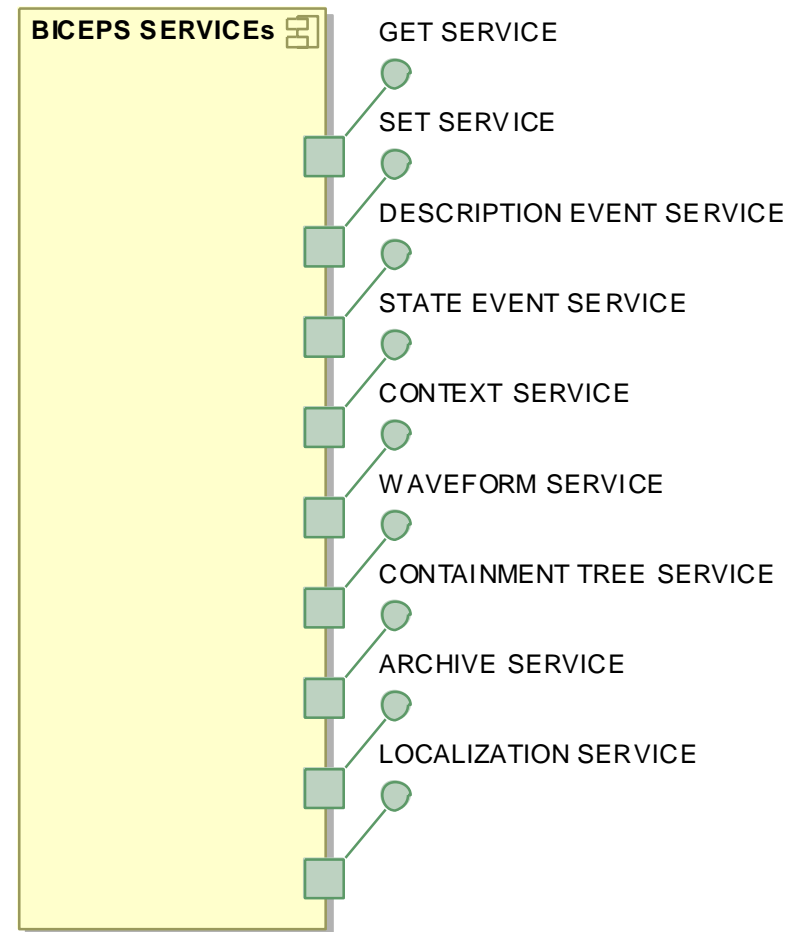

Messages

Separation of operations and messages


- A service operation can be considered as an interprocess function taking input parameters and computing output results
- The BICEPS Message Model defines the input parameters and output results of the BICEPS Service Model
- Naming convention for messages:
 - Request: <BaseName>
 - Response: <BaseName> + Response
 - Notification: <BaseName> + Report

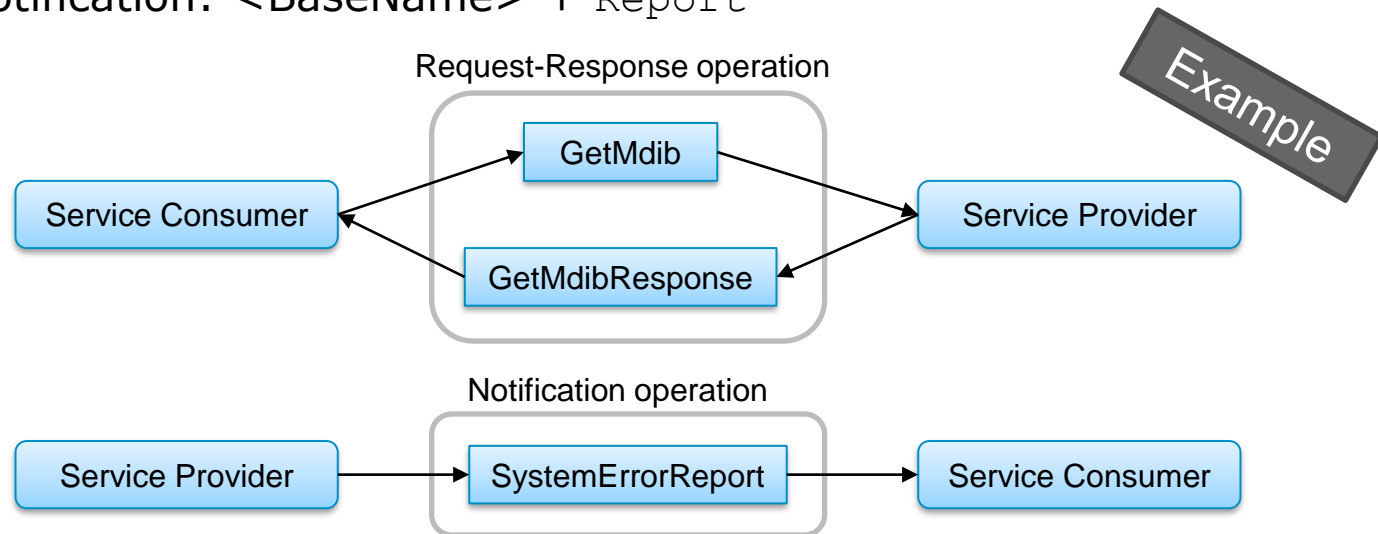


Messages

Two report flavors

- Episodic reports: delivered on change
- Periodic reports: delivered continuously at a given period
- Naming convention for reports that are available either episodically or periodically:
 - `[Periodic | Episodic] + <BaseName> + Report`

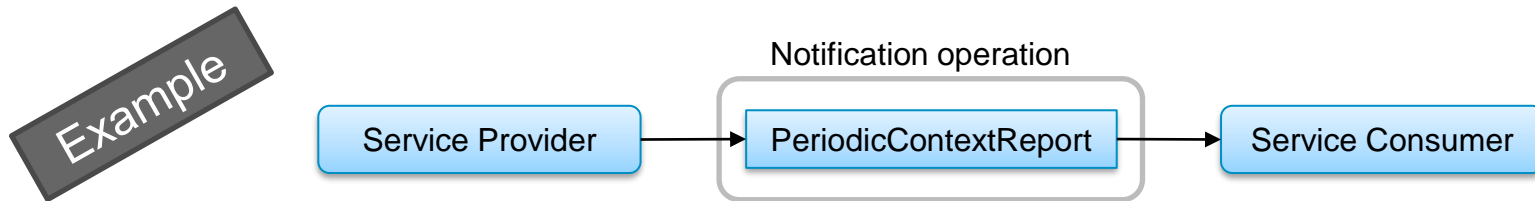


Component view

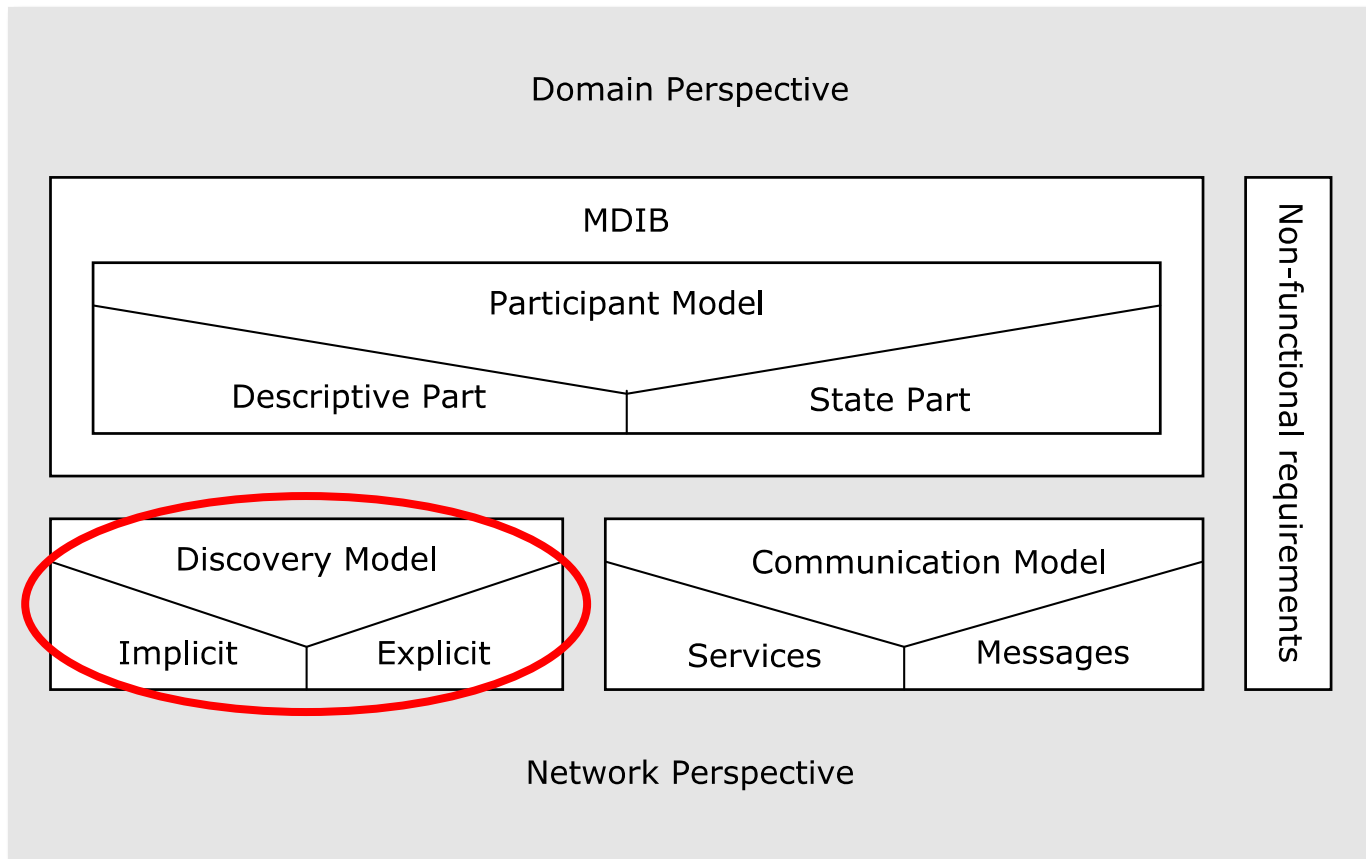


Discovery model

Explicit discovery

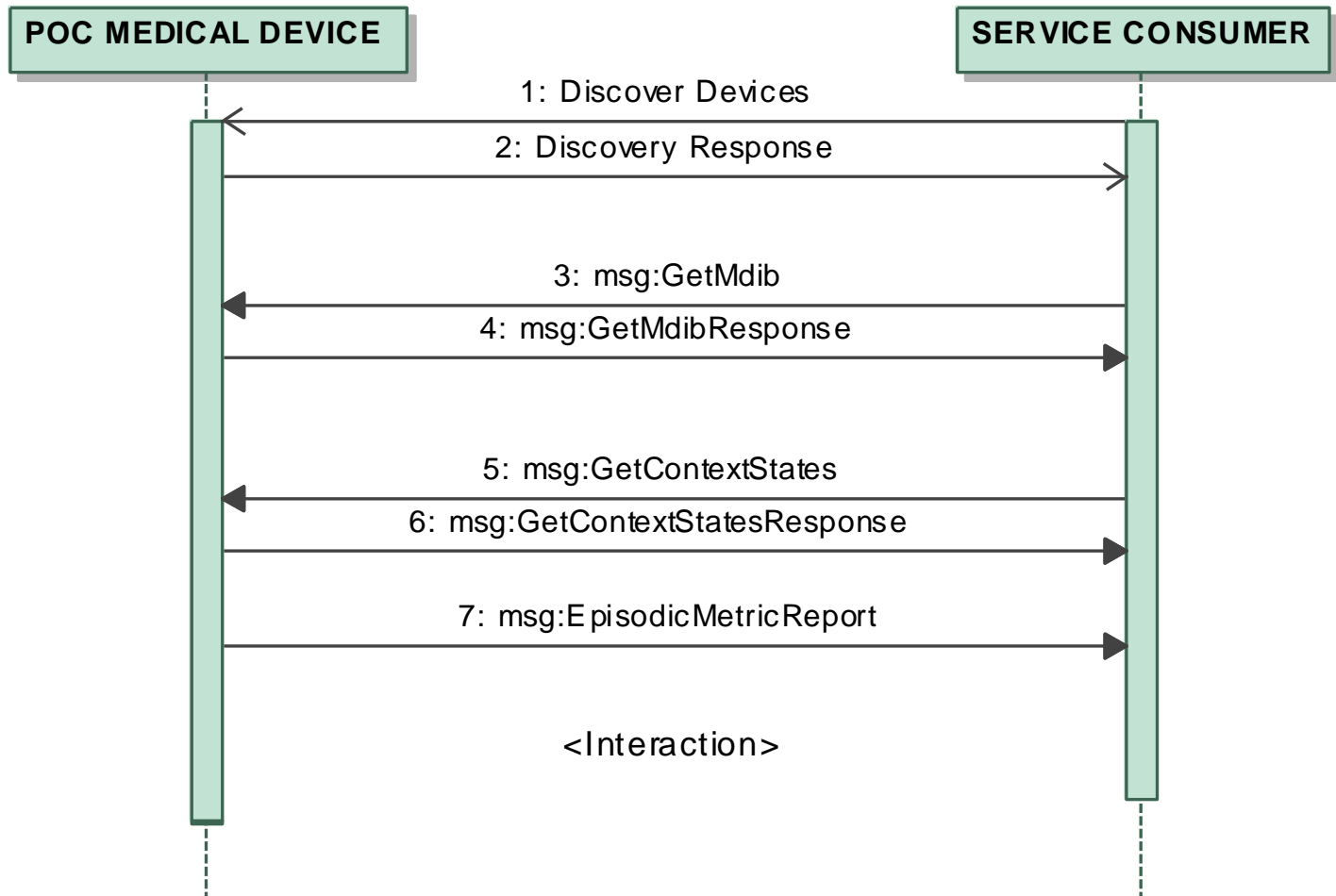


Discovery model

Implicit discovery

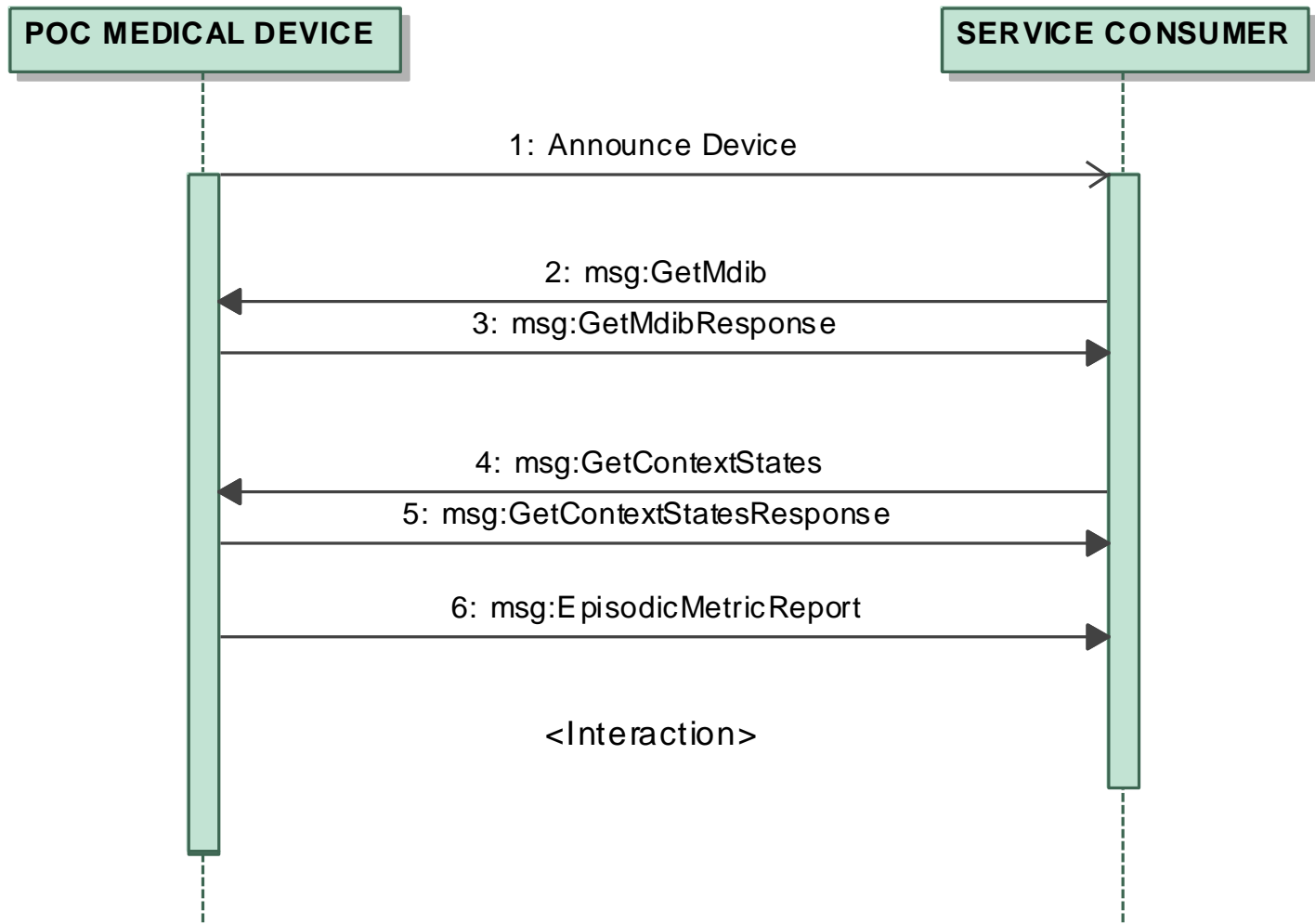

Component view

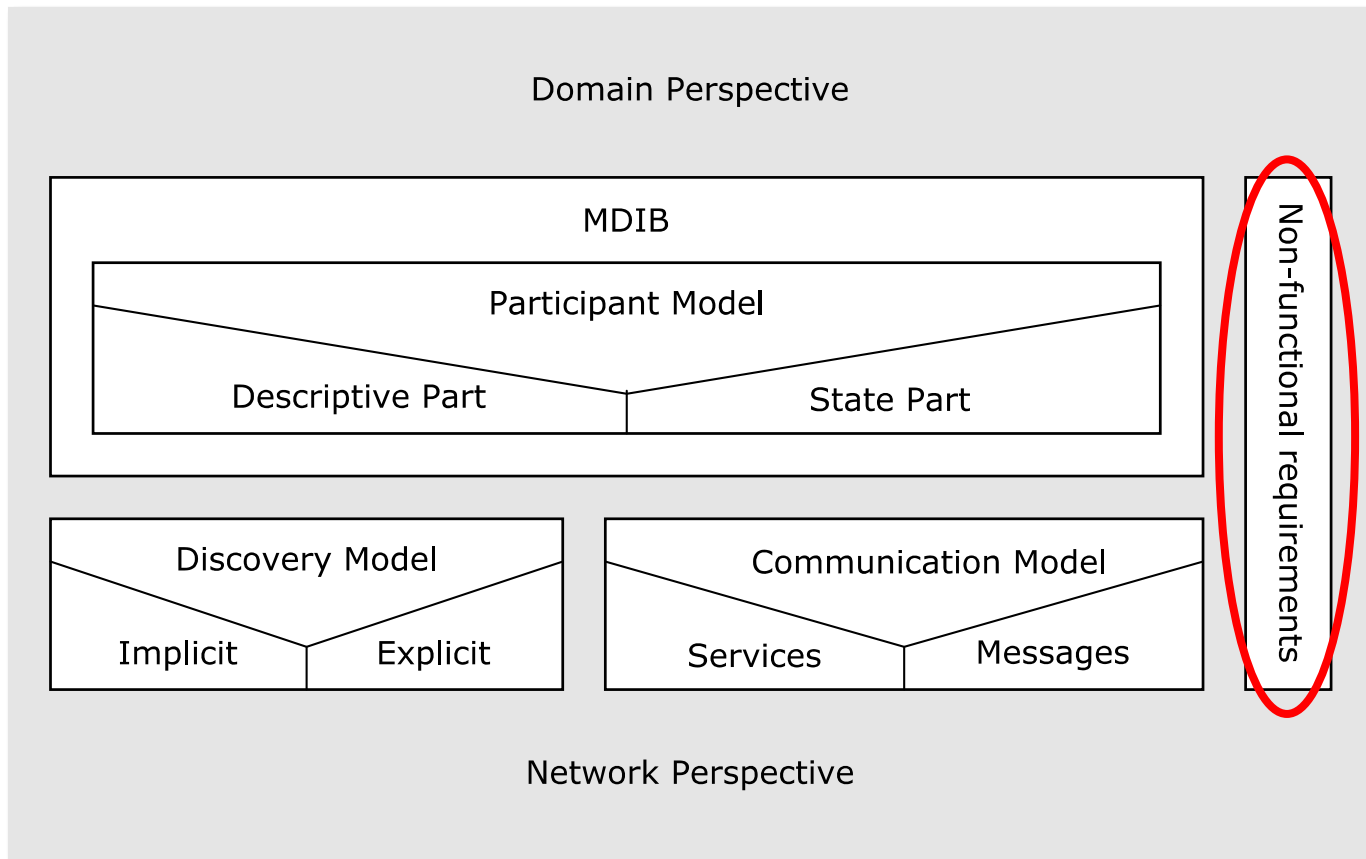

Non-functional Requirements

- Patient Safety & Cyber Security Concerns
 - A BICEPS compliant binding is required to be able to establish
 - confidentiality between participants
 - trust through authorization
 - data integrity
 - connection loss detection between participants that exchange messages
 - accountability between participants that exchange messages
- Clinical Effectiveness & Regulatory Considerations
 - A BICEPS compliant binding is required to
 - support time synchronization between participants
 - support provision of defining QoS metrics between participants
 - distinguish unique messages in a sequence of messages with potential duplicates

Thank you for your attention!

Contact information

David Gregorczyk

david.gregorczyk@draeger.com